


Facein™ Combo



SANMED
ACTIVE HYGIENE

Surface, Equipment & Environment Disinfectant Concentrate

EFFECTIVE AGAINST HIV & HBV.

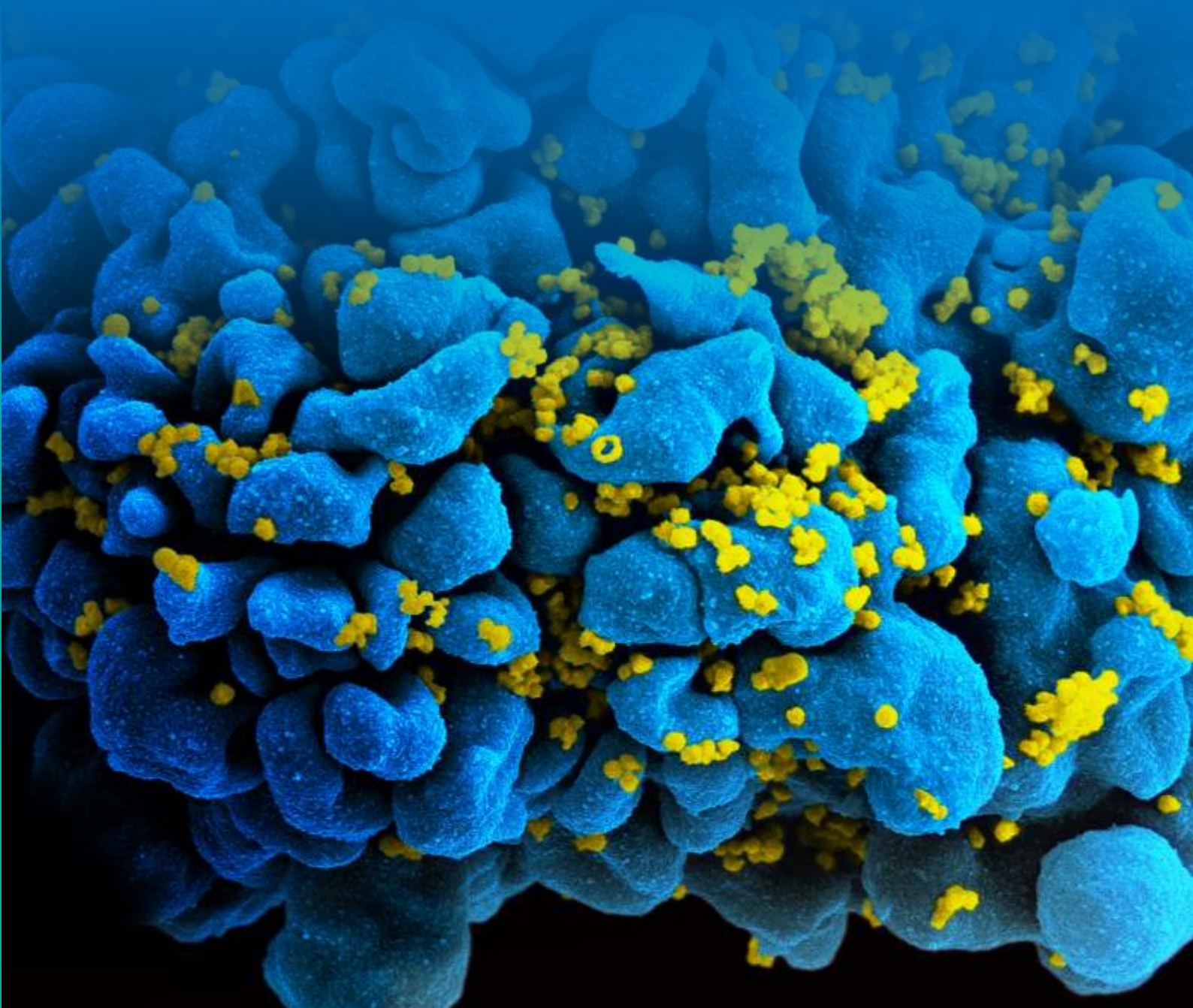

Facein™ Combo

SURFACE, EQUIPMENT & ENVIRONMENT DISINFECTANT CONCENTRATE

EFFECTIVE AGAINST HIV & HBV.

Composition

Each 100 g contains.

1, 6 Dihydroxy, 2, 5-Dioxahexane (Chemically bound Formaldehyde)	- 11.2 g
Glutaraldehyde USP	- 5.0 g
Benzalkonium Chloride Soln IP (Eq. To Benzalkonium Chloride)	- 5.0 g
Alkyl Urea Derivative	- 3.0 g

Salient Features

- ❖ Excellent Virucidal activity against HIV, HBV, Polio Virus etc.
- ❖ Broad spectrum Bactericidal (including MRSA), Sporocidal, Tuberculocidal and Fungicidal
- ❖ Corrosion inhibitors for greater safety
- ❖ Sustained high degree of asepsis.

FACEIN™ COMBO Ensures

- ❖ Compatibility with all types of surfaces
- ❖ Does not irritate eyes
- ❖ Economical – 500 ml =100 Liters prepared solution (0.5%)
- ❖ Excellent Cleansing
- ❖ Maintenance of lowest microbial count
- ❖ No corrosion on frequent use
- ❖ Pleasant odour
- ❖ Prolonged Residual action
- ❖ Rapid disinfection within 10-30 minutes.

Instructions for Use

- ❖ **High-risk Areas:** ICU, ICCU, Oncology dept., OTs, Minor OTs, Pharmaceutical Sterile sections, Neo-Natal dept., Trauma Care, AKD etc. Spray over all exposed surfaces with a hand compression (pneumatic) sprayer
- ❖ **Non-critical Areas:** For Emergency disinfections and between cases in wards, corridors, OPDs, Clinics, Mobile Units, Laboratories, Dental Clinics etc.
- ❖ Add 25 ml to 50 ml of FACEIN COMBO to 10 Ltr. (2/3 bucket) of tap water to get 0.25% to 0.5 % solution.



SANMED
ACTIVE HYGIENE

Sanmed Healthcare Pvt Ltd.,

401, Abhis Hiranya Guttala Begumpet Madhapur Hyderabad 500081 India
Telephone : +91 40 33574705
Email ID : info@sanmed.in

Factory: Plot No: 56, Biotech Park-3, Karakapatta, Medak Dist, A.P, India



- ❖ Spray or mop liberally on exposed surfaces, keeping wet for 10-20 mins for optimum result.
- ❖ FACEIN™ COMBO can be used in OT Care Machine when used in 1 % or less than 1% dilution.

Precautions

- ❖ Turn off air Conditioners and Fans.
Do not spray without wearing mask.
- ❖ Avoid contact with concentrated solution
- ❖ Keep out of reach of Childers
- ❖ Do not mix with detergents or other disinfectants.

Standard Guidelines for Disinfection

- ❖ **At 2% weekly:** Add 100ml FACEIN™ COMBO in 5 Ltr of water
- ❖ **At 1% alternate days:** Add 50ml FACEIN™ COMBO in 5Ltr of water
- ❖ **At 0.5% Daily:** Add 25ml FACEIN™ COMBO in 5 Ltr of water

Note

- ❖ To reduce evaporation before and after application, close door and windows and switch off fans and Ac's.
- ❖ Disinfected surfaces must remain in contact i.e wet with solution for 10 to 20 minutes for optimum result.

Packing

- ❖ 500 ml. ❖ 5 Ltr.

Mopping is recommended over Fogging (CDC Guidelines - U.S.A. Canadian Hospital Infection Control Guidelines). If using, Fogging machines, use 1 % FACEIN™ COMBO daily i.e. 20 ml in 2 lts of water.

www.sanmed.in